





Closed IV Catheter Device


WinFusion

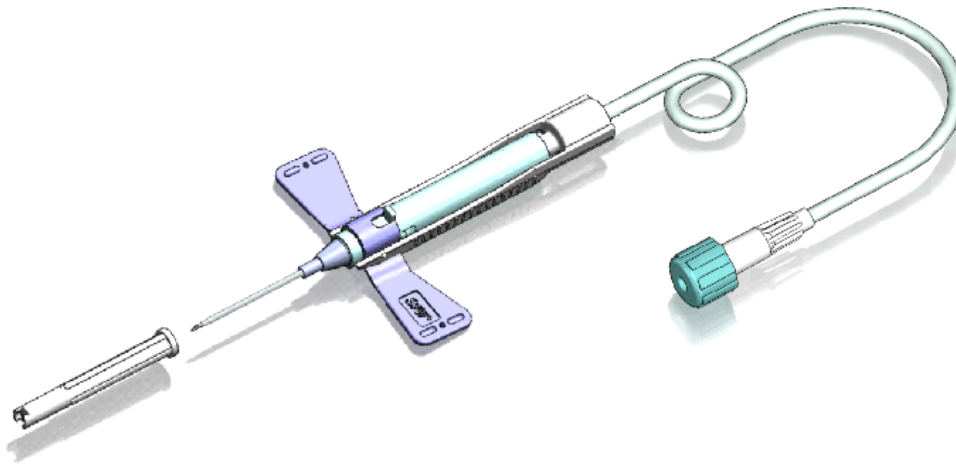

Closed IV catheter with the very unique design offering solutions against the existing risks

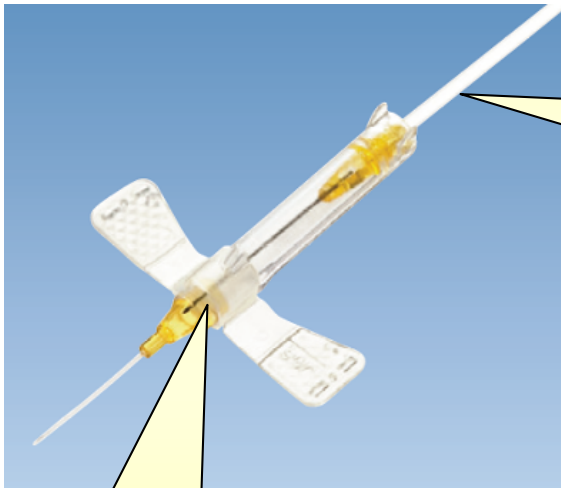

Tubing is integrated for creating the closed system

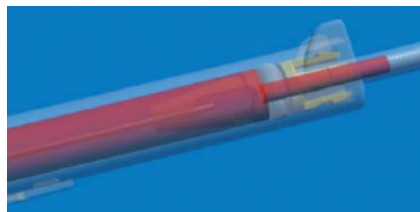


The closed system prevents blood exposure and hub contamination.

Needle never comes out

The needle is retained inside the catheter body.

- No risk of Needle-stick injury
- No need to dispose of the needle separately



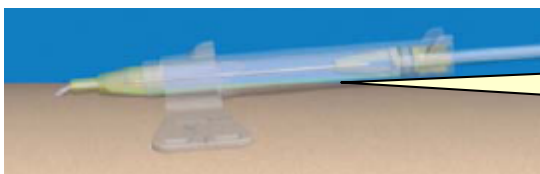


Flat and Seamless shape


- Gives gentle fixation on the skin.
- No risk of drug exposure.

Operation Procedure


Three ways of Holding Position of WinFusion

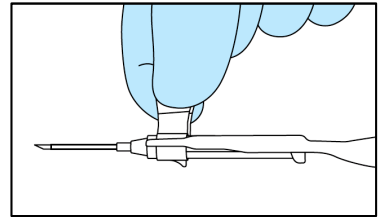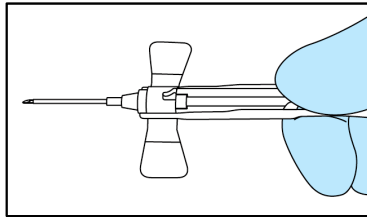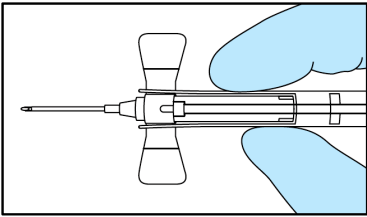


1. Venipuncture

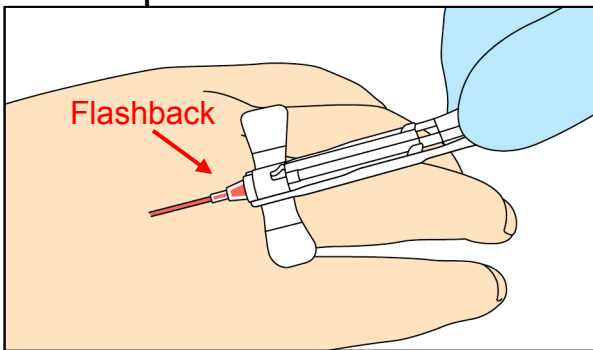

2. Indwelling

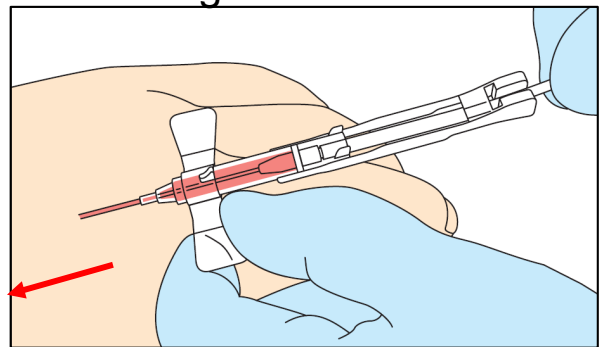

3. Retraction

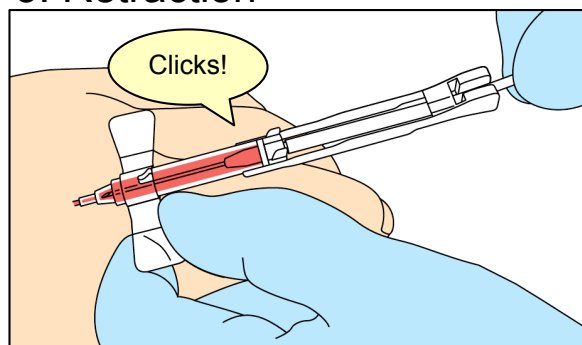

Product Configuration

| Gauge | Length | Hub Color | Packaging |
|-------|-------------|-----------|--------------------|
| 20G | 3/4" (19mm) | Pink | 50 x 20 = 1 carton |
| 22G | 3/4" (19mm) | Blue | 50 x 20 = 1 carton |
| 24G | 3/4" (19mm) | Yellow | 50 x 20 = 1 carton |



JMS CO.,LTD.
 12-17 Kako-machi, Naka-ku
 Hiroshima, 730-8652 Japan